
SCULPTURAL SCULPTURAL BRUSH

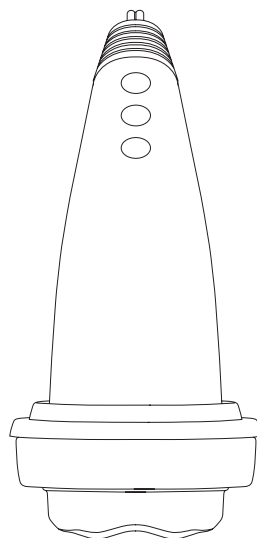- ES.** Consejos y advertencias de seguridad
EN. Safety advice and warnings
FR. Conseils et avertissements de sécurité
PT. Conselhos e advertências de segurança
IT. Consigli e avvertimenti di sicurezza
CA. Consells i advertències de seguretat
DE. Ratschläge und sicherheitshinweise
NL. Veiligheidsadvies en waarschuwingen
RO. Recomandări siguranță și avertizări
PL. Wskazówki i ostrzeżenia bezpieczeństwa
BG. Съвети и предупреждения за безопасност
EL. Συμβουλές και προειδοποιήσεις ασφαλείας
RU. Советы по безопасности и предупреждения
DA. Sikkerhedsråd og advarsler
NO. Sikkerhetsråd og advarsler
SV. Säkerhetsråd och varningar
FI. Urvallisuusohejeet ja varoitukset
TR. Güvenlik tavsiyeleri ve uyarılar
HE. עצות ואזהרות בטיחות
AR. نصائح وتحذيرات تتعلق بالسلامة

solac

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**MASSAGER
SCULPTURAL
SCULPTURAL BRUSH**

Dear customer,

Many thanks for choosing to purchase a SOLAC brand product.

Thanks to its technology, design and operation and the fact that it exceeds the strictest quality standards, a fully satisfactory use and long product life can be assured.

SAFETY ADVICE AND WARNINGS

- Read these instructions carefully before switching on the appliance and keep them for future reference. Failure to follow and observe these instructions could lead to an accident.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- This appliance is for household use only, not professional or industrial use.

- Ensure that the voltage indicated on the rating label matches the mains voltage before plugging in the appliance.
- Connect the appliance to a mains socket rated at least 10 amperes.
- The appliance's plug must fit into the mains socket properly. Do not alter the plug. Do not use plug adaptors.
- Do not use the appliance if the supply cord or the plug is damaged.
- If the appliance enclosure breaks, immediately disconnect the appliance from the mains to prevent the possibility of an electric shock.
- Do not pull on the supply cord. Never use the supply cord to lift up, carry or unplug the appliance.
- Do not wrap the supply cord around the appliance.
- Check the state of the supply cord. If it is damaged it will increase the risk of electric shock.
- Do not allow the supply cord to hang freely or to come into contact with the appliance's hot surfaces.
- Do not touch the plug with wet hands.
- Do not use or store the appliance outdoors.
- Do not leave the appliance out in the rain or exposed to moisture. If water gets into the appliance, this will increase the risk of electric shock.
- The temperature of the accessible surfaces may be high when the appliance is in use.
- Do not touch any heated part, as this may cause serious burns.

SERVICE

- Any misuse or failure to follow the instructions for use renders the guarantee and the manufacturer's liability null and void.

ANOMALIES AND REPAIR

- Take the appliance to an authorised technical support service if problems arise. Do not try to dismantle or repair without assistance, as this may be dangerous.

**FOR EU PRODUCT VERSIONS AND/OR IN THE CASE
THAT IT IS REQUESTED IN YOUR COUNTRY:
ECOLOGY AND RECYCLABILITY OF THE PRODUCT**

- The materials of which the packaging of this appliance consists are included in a collection, classification and recycling system. Should you wish to dispose of them, use the appropriate public containers for each type of material.
- The product does not contain concentrations of substances that could be considered harmful to the environment.

This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an authorised waste agent for the selective collection of waste electrical and electronic equipment

(WEEE).

- This appliance complies with Directive 2014/35/EU on Low Voltage, Directive 2014/30/EU on Electromagnetic Compatibility, Directive 2011/65/EU on the restrictions of the use of certain hazardous substances in electrical and electronic equipment and Directive 2009/125/EC on the ecodesign requirements for energy-related products.

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ENGLISH

WARRANTY AND TECHNICAL ASSISTANCE

- This product enjoys the recognition and protection of the legal guarantee in accordance with current legislation. To enforce your rights or interests you must go to any of our official technical assistance services.
- You can find the closest one by accessing the following web link: <http://solac.com/>
- You can also request related information by contacting us.
- You can download this instruction manual and its updates at <http://solac.com>

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MASSAGER SCULPTURAL SCULPTURAL BRUSH

DESCRIPTION

- A Main body
- B ON/OFF and heat button
- C Speed selector
- D Mode selector
- E Supply cord
- F Red LEDs
- G Massage head adaptor
- H Wave shaped head
- I Finger massage head
- J Serrated roller head
- K Exfoliating head **
- L Smooth roller head
- M Ball massage head *
- N Parallel spiky disk massage head *
- O Fabric cover
- P Storage bag

(*) Only available in the model Sculptural

(**) Only available in the mode Sculptural Brush

If the model of your appliance does not have the accessories described above, they can also be bought separately from the Technical Assistance Service.

USE AND CARE

- Do not use the massager on the abdomen area if you are pregnant or believe that you might be pregnant.
- Should you suffer from any health problems, e.g., pace maker fitted, consult your doctor.
- Do not use the massager directly on:
 - Skin rashes or inflammation.
 - Burns, open wounds, infected skin.
- Do not operate the device too close to the hair, tie up the hair before use to prevent the hair being entangled.
- Do not use it if you suffer from any of the following illnesses or ailments: varicose veins, poor circulation, bruising, and inflammation of the veins, thrombosis, torn skin or diabetes.
- Do not use the appliance if you are taking medication that may reduce sensory capacity.

- Do not touch any of the appliance's moving parts while it is operating.
- Fully extend the appliance's supply cord before each use.
- Do not use the appliance if its accessories are not properly fitted.
- Do not use the appliance if the accessories attached to it are defective. Replace them immediately.
- Do not use the appliance if the on/off button does not work.
- If you use the appliance in a bathroom or similar place, unplug the appliance from the mains when it is not in use, even if it is only for a short while, as proximity to water poses a risk, even if the appliance is switched off.
- Disconnect the appliance from the mains when not in use, before cleaning, making any adjustment or accessory change.
- This appliance is for household use only, not professional or industrial use.
- Store this appliance out of the reach of children and/or those with reduced physical, sensory or mental abilities or those unfamiliar with its use.
- Do not store the appliance if it is still hot.
- Keep the appliance in good condition. Check that the moving parts are not misaligned or locked and make sure there are no broken parts or anomalies that may prevent the appliance from operating correctly.
- Use the appliance and its accessories and tools in accordance with these instructions, taking into account the working conditions and the work to be performed. Use the appliance for operations different from those intended could result in a hazardous situation.
- Never leave the appliance connected and unattended if it is not in use. This saves energy and prolongs the life of the appliance.
- Never rest the appliance on a surface while in use.
- Do not use the appliance for more than 20 minutes at any one time. Allow to cool for one hour before subsequent use.
- Do not use the appliance on pets or animals.
- CAUTION: Do not fall asleep while using this appliance, as this could be harmful.

INSTRUCTIONS FOR USE

BEFORE USE

- Make sure that all product's packaging has been removed.
- Please read the "Safety advice and warnings" booklet carefully before first use.

GENERAL RECOMMENDATIONS FOR USE

- The intensity control allows you to regulate the intensity of the massage as required. Start with a gentle speed and progressively increase depending on your sensitivity to the massage.
- The red light is an additional massage function. Use it to strengthen the effectiveness of the massage.
- If you use specific treatment creams apply these after the massage. Prepare the skin and tissues for the absorption of the ingredients. The cream can be applied using the "Wave shaped head" (H) or the "Smooth roller head" (L)
- If you use the "Wave shaped head" (H) head protection attachment (O), do not apply creams.

RECOMMENDED DURATION OF TREATMENT

- The treatment can last between 15-20 minutes.
- The appliance will automatically stop after 20 minutes, if you wish to continue using the appliance after this time, allow it to cool down before continuing.
- For a more effective treatment, use the appliance regularly. Use the appliance each day until you notice the results, you can then use it to maintain these results. To maintain the results over the long-term, adapt your lifestyle habits and food.

USE

- Remove the protective cover (O).
- Remove the default massage head attached. (Fig. 1)
- Choose a massage head (see "Choosing a massage head").
- You may alternate the use of different heads during the same massage.
- Fit the head corresponding to the desired massage onto the head adaptor. (Fig. 2)
- Extend the supply cord completely before plugging it in.
- Connect the appliance to the mains.
- The led light of the ON/OFF button (B) will flash, indicating that the appliance is in stand-by mode.
- Hold the appliance with the control panel facing you, holding it by the anti-slip area in the back (Fig. 3).
- Turn the appliance on, by using the ON/OFF button (B).
- Press and hold the ON/OFF button (B) to activate the heating function.
- To avoid a fierce start, the appliance is turned on in normal operation and super low speed.
- Select the desired massage mode, by using the speed and mode selector (C) (D) (see "Choosing a pre-set mode") and place the massager on the area on which you are going to work.

- Hold the appliance firmly with both hands and move it over the skin in the direction of the blood circulation (see "Massage Direction").
- Without pressing too hard, simply glide the massager over the area to be treated. Pressing the appliance too hard against the skin will prevent the motor from operating efficiently.

CHOOSING A MESSAGE HEAD

- FINGER-PRINT (I): Finger massage head. Moving pressure tips.
- MAX-FORCE (J, N): Serrated roller heads. Moving, vibrating helixes. This tones the muscles, reduces volume, warms up the muscles, relaxes and relieves stress, for example.
- PRO-WAVE (H): Wave shaped head. Tones and improves blood circulation. For use on the thigh and buttock areas.
- ROLL-MASSAGE (L): Smooth roller head. Sequence of long, sliding movements. Relaxes and improves blood circulation.
- DRY BRUSH (K): Exfoliating brush, removes dead skin and stimulates blood circulation.
- BALL-MASSAGE (M): Ball massage head. Stimulates and tones the muscles, helping to firm the skin.
- FABRIC COVER (O): For exclusive use with the PRO-WAVE head (H). For greater comfort, cover the head with the fabric cover. Never use it with creams. (Fig. 4)

CHOOSING A PRESET MODE

- NORMAL OPERATION with 4 speeds.
- Select the desired speed by pressing the speed selector button (C):
 - Super low speed: ~1.800 RPM
 - Low speed: ~2.200 RPM
 - Medium speed: ~2.600 RPM
 - High speed: ~3.000 RPM

FIXED PRESET MODES

- Select the desired fixed pre-set mode by pressing the mode selector button (D):
 - STIMULATING MODE: 3 automatic frequency variation per minute (120, 150 and 200 beats).
 - SYSTEMATIC WAVE MODE: cycle time of 8 seconds between super low speed and high speed.
 - FIRING MODE: fixed 120 beats per minute.
 - FAST FIRING MODE: fixed 200 beats per minute.
- To go back to normal operation, press the speed selector button (C).

MESSAGE DIRECTION

- ALWAYS glide the appliance over the skin in the direction of the blood circulation. In an upward direction.
- NEVER directly massage the spine.
- NEVER directly massage a bone (knee, ankle, and elbow). The protruding parts of the massage head could hit the bone, causing an injury.
 - Legs: Place the appliance on the thigh, close to the knee, and glide it gently upwards towards the hip. (Fig. 5).
 - Buttocks: Combine upward and circular movements. (Fig. 6).
 - Abdomen: Place the appliance on the navel and move it gently outwards. (Fig. 7). Never make circular movements around the navel.
 - Arms: Place the appliance close to the elbow and glide it gently towards the armpit. Repeat the massage from the armpit towards the elbow. (Fig. 8).
 - Back: Place the appliance in the area of the spine and glide it gently upwards towards the shoulder blade. (Fig. 9).

SKIN PREPARATION – EXFOLIATING

- It is essential to start by preparing the skin properly before performing the massage in the area where it is required, so we should exfoliate the skin using the natural bristle brush (DRY BRUSH) (K), which will leave the area free of dead cells and make the skin softer. Avoid damaging the skin through excessive exfoliation (once per week is sufficient). Gently massage the entire area for 5 minutes. This step helps to stimulate your blood circulation and metabolism.
- The skin should always be brushed toward the heart, starting with the feet and moving up to the hips in the case of legs, and from wrist to shoulder for arms. Persist on rougher areas and be careful with sensitive areas. Adjust the intensity of exfoliation based on the condition of the skin, the body area and particular sensitivity.

ONCE YOU HAVE FINISHED USING THE APPLIANCE

- Turn the appliance off using the ON/OFF button (B).
- Unplug the appliance to the mains.
- Clean the appliance.

SAFETY THERMAL PROTECTOR

- The appliance has a safety device, which protects the appliance from overheating.
- If the appliance turns itself off, leave the appliance to cool down before attempting to restart it. If you wish to continue using the appliance, allow it to cool down before proceeding.

CLEANING

- Disconnect the appliance from the mains and let it cool before undertaking any cleaning task.
- Clean the equipment using a damp cloth with a few drops of washing-up liquid and then dry.
- Do not use solvents, or products with an acid or base pH such as bleach, or abrasive products, for cleaning the appliance.
- Do not immerse the appliance in water or any other liquid, or place it under a running tap.

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