
AQUA PULSE

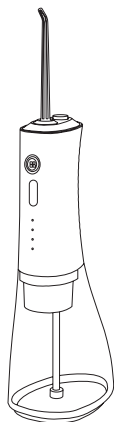- ES.** Consejos y advertencias de seguridad
- EN.** Safety advice and warnings
- FR.** Conseils et avertissements de sécurité
- PT.** Conselhos e advertências de segurança
- IT.** Consigli e avvertimenti di sicurezza
- CA.** Consells i advertències de seguretat
- DE.** Ratschläge und sicherheitshinweise
- RO.** Recomandări siguranță și avertizări
- PL.** Wskazówki i ostrzeżenia bezpieczeństwa
- BG.** Съвети и предупреждения за безопасност
- AR.** نصائح وتحذيرات تتعلق بالسلامة

solac

ES

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**DENTAL IRRIGATOR
AQUA PULSE**

Dear customer,

Many thanks for choosing to purchase a SOLAC brand product.

Thanks to its technology, design and operation and the fact that it exceeds the strictest quality standards, a fully satisfactory use and long product life can be assured.

**SAFETY ADVICE AND
WARNINGS**

Read these
instructions

carefully before
switching on the appliance
and keep them for future
reference. Failure to follow
and observe these
instructions could lead to
an accident.

- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory

or mental capabilities or lack of experience and knowledge under supervision.

- This oral irrigator shall be used only with cold or warm water (below 40°C).
- Children shall not play with the appliance.
- This appliance is for household use only, not professional or industrial use.

This appliance
has a detachable
supply unit.
(ZD006M050100EU).

- This appliance is only to be used with the power supply unit provided with the appliance.
- This appliance contains batteries that are non-replaceable. When the

battery is at end of life, the appliance shall be properly disposed of.

- **WARNING:** For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
- Ensure that the voltage indicated on the rating label matches the mains voltage before plugging in the appliance.
- The appliance's plug must fit into the mains socket properly. Do not alter the plug. Do not use plug adaptors.
- Do not use the appliance if the detachable supply unit is damaged.
- If the appliance enclosure breaks, immediately disconnect the appliance from the mains to prevent the possibility of an electric shock.
- Do not touch the plug with wet hands.

SERVICE

- Any misuse or failure to follow the instructions for use renders the guarantee and the manufacturer's liability null and void.

ANOMALIES AND REPAIR

- Take the appliance to an authorised technical support service if problems arise. Do not try to dismantle or repair

without assistance, as this may be dangerous.

FOR EU PRODUCT VERSIONS AND/OR IN THE CASE THAT IT IS REQUESTED IN YOUR COUNTRY: ECOLOGY AND RECYCLABILITY OF THE PRODUCT

- The materials of which the packaging of this appliance consists are included in a collection, classification and recycling system. Should you wish to dispose of them, use the appropriate public containers for each type of material.
- The product does not contain concentrations of substances that could be considered harmful to the environment.

This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an authorised waste agent for the selective collection of waste electrical and electronic equipment (WEEE).

- This appliance complies with Directive 2014/35/EU on Low Voltage, Directive 2014/30/EU on Electromagnetic Compatibility, Directive 2011/65/EU on the restrictions of the use of certain hazardous substances in electrical and electronic equipment and Directive 2009/125/EC on the ecodesign requirements for energy-related products.
- The following information is regarding regulation 2019/1782 on ecodesign

requirements for external power
supplies:

Trademark and address	Solac Small Electro, S.L. Avda. Barcelona, s/n, 25790 Oliana, Lleida (Spain)
Model	ZD006M050100EU
Input voltage	100-240 V
Input AC frequency	50/60 Hz
Output voltage	5,0 V
Output current	1,0 A
Output power	5,0 W
Average active efficiency	74,00%
Efficiency at low load (10 %)	N.A.
No-load power consumption	0,09 W

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ENGLISH

WARRANTY AND TECHNICAL ASSISTANCE

- This product enjoys the recognition and protection of the legal guarantee in accordance with current legislation. To enforce your rights or interests you must go to any of our official technical assistance services.
- You can find the closest one by accessing the following web link: <http://solac.com/>
- You can also request related information by contacting us.

- You can download this instruction manual and its updates at <http://solac.com/>

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English (Original instructions)

**DENTAL IRRIGATOR
AQUA PULSE****DESCRIPTION**

- A Nozzle
- B Nozzle release button
- C ON/OFF button
- D Mode button
- E Mode pilot lights
- F Battery pilot light
- G Removable water tank
- H Water inlet lid
- I Charger
- J Suction tube
- K Filter
- L Charge socket

If the model of your appliance does not have the accessories described above, they can also be bought separately from the Technical Assistance Service.

USE AND CARE

- Do not use the appliance if the parts or accessories are not properly fitted.
- Do not use the appliance if the accessories or consumables attached to it are defective. Replace them immediately.
- Do not switch the appliance on without water.
- Do not use the appliance if the on/off button does not work.
- Do not use the appliance if it is tipped up and do not turn it over.
- Do not turn the appliance over while it is in use or connected to the charger.
- Respect the MAX level indicator.
- If you use the appliance in a bathroom or similar place, unplug the appliance from the mains when it is not in use, even if it is only for a short while, as proximity to water involves a risk, even if the appliance is disconnected.
- Unplug the appliance from the mains when not in use and before undertaking any cleaning task, reclape and accessory or adjustment.
- Unplug the appliance from the mains before refilling the water tank.
- Keep the appliance in a dry, dust-free place, out of direct sunlight.
- Keep the appliance in good condition. Check that the moving parts are not misaligned or jammed and make sure there are no broken parts or anomalies that may prevent the appliance from operating correctly.
- Use the appliance and its accessories and tools in accordance with these instructions, taking into account the working conditions and the work to be performed. Use the appliance for operations different from those intended could result in a hazardous situation.
- Do not open the battery/ies under any circumstances whatsoever.

- Never leave the appliance connected and unattended if it is not in use. This saves energy and prolongs the life of the appliance.
- Never rest the appliance into a surface while in use.
- The appliance has a 2 minutes timer function, which will automatically shut down after continuous use for 2 minutes. If needs to continue using it, simply press the on /off button again.
- When the water tank is fully filled with 320ml of water, the duration of use may vary depending on the cleaning mode.
- Do not use the appliance with pets or animals.
- Recharge the battery/ies only with the charger supplied by the manufacturer. A charger that is suitable for one type of battery may cause explosion or fire hazard when used with other battery/ies.
- Use the device only with the battery/ies for which has been specifically designed. The use of any other battery/ies can cause an explosion or fire hazard.

CHARGING THE BATTERY

Caution: A new battery is not fully charged and must therefore be charged before the appliance is used for the first time.

- Charge the battery at an atmospheric temperature of between 0 and 45°C.
- Check and ensure that the product and hands are in a dry state.
- Plug the charger into the charging port to the mains.
- A pilot light will come on to indicate that the charger is ready to begin charging the battery.
- Approximately 6 hours are needed for a complete charge of the battery. Never attempt to over charge the battery.
- Connect the charger cable to the jack using the socket on the appliance. The socket is located at the back part of the appliance.
- The pilot light will come on and flicker to indicate that the charger cable is charging the battery.
- At the end of the time period specified above, remove the charger cable from the appliance.
- When the battery is charged disconnect the charger from the mains supply.

NOTE: The appliance can't be used while it is connected to the charger.

HOW TO CHECK THE BATTERY'S CHARGE LEVEL

- If pilot light flicker this means that the battery has very little power remaining.
- When the appliance is connected to the mains, the pilot light will flash orange color if the battery level is lower than 20% or green if it's between 20 -100%.
- The pilot light lights up for 30 seconds and then turns off, when battery is fully charged.
- In use, the pilot lights in blue if the battery level is between 5-100% or flashes in orange if the battery level is lower than 5%.
- When the battery level reaches 0% the pilot light flashes quickly in orange.

NOTE: For a long battery life, it is advisable to charge the battery, only, when it has run out completely.

INSTRUCTIONS FOR USE BEFORE USE

- Make sure that all the product's packaging has been removed.
- Please read carefully the "Safety advice and warnings" booklet before first use
- Before using this product for the first time, it is important to make a full charge of battery.

FILLING WITH WATER

- You must fill the water tank before turning on the appliance.
- Hold the appliance horizontally.
- Open the water inlet lid and fill the water tank with an appropriate amount of water taking care to observe the MAX level.
- Close the water inlet lid.
- It is recommended to fill the water tank with warm water below 40°C or running tap water.
- This appliance is designed to operate with water only, do not fill it with any other liquid or substance.

FITTING THE NOZZLE

- The appliance is designed to be used with an interchangeable nozzle fitted in.
- To attach a nozzle to it, simply press the nozzle in its housing.
- To remove it, press the nozzle release button and pull the nozzle.

Note: Do not press the nozzle release button during use to avoid accidents or injuries.

USE

- Don't use the appliance if it is connected to the main.
- Hold the appliance vertically, lean over the sink and insert the nozzle into the mouth.

- Guide the appliance to direct the flow of water in the desired direction.
- Turn the appliance on, by short pressing the on/off button.
- The pilot light comes on and stay on.
- Select the desired working mode.
- The appliance stops working automatically after 2 minutes, but you can stop it manually at any time by pressing the on/off button again.

ONCE YOU HAVE FINISHED USING THE APPLIANCE

- Turn the appliance off, using the on/off button.
- Remove the water from the inside of the appliance.
- Clean the appliance.

CLEANING

- Disconnect the appliance from the mains and allow it to cool before undertaking any cleaning task.
- Clean the equipment with a damp cloth with a few drops of washing-up liquid and then dry.
- Do not use solvents, or products with an acid or base pH such as bleach, or abrasive products, for cleaning the appliance.

HOW TO REMOVE THE BATTERIES FROM THE INTERIOR OF THE DEVICE

- This appliance contains batteries that are only replaceable by skilled persons.
- In case you need to replace the batteries, take the appliance to an authorised technical support service. Do not try to dismantle or repair without assistance, as this may be dangerous.

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